

Topbar

Background

CAT Prompts are created in CAT4 (on the practice level) or in PAT CAT (on the PHN or other organisation level). Prompts created in PAT CAT will have to be activated by the practice before they are displayed in Topbar and are sent once a day (over night) from PAT CAT to Topbar. Prompts created in CAT4 in the practice are effective immediately and don't need any further activation in Topbar to be displayed.

Topbar has to be linked to CAT4 in the practice to display CAT Prompts created locally, the details are explained [here](#).


The clinic has to be set up in PAT CAT to receive prompts created in PAT CAT, the details are explained in the [PAT CAT Prompts guide](#).

Displaying Prompts

Once a patient is opened in the clinical system matching the criteria of the prompt, the CAT Prompts App in Topbar will display a notification. As with all Topbar apps, the notification is displayed as a red circle with the number of notifications in the circle. In the example below the patient has two notification in the CAT Prompts App.



To open the CAT Prompts App click on 'cat plus prompts' in Topbar to display the full screen. Prompts displayed for the patient are labeled according to the name given to the prompt when created, and also show which software was used to create the prompt. In the example below we can see the same prompt originating once from PAT CAT and once from CAT4.

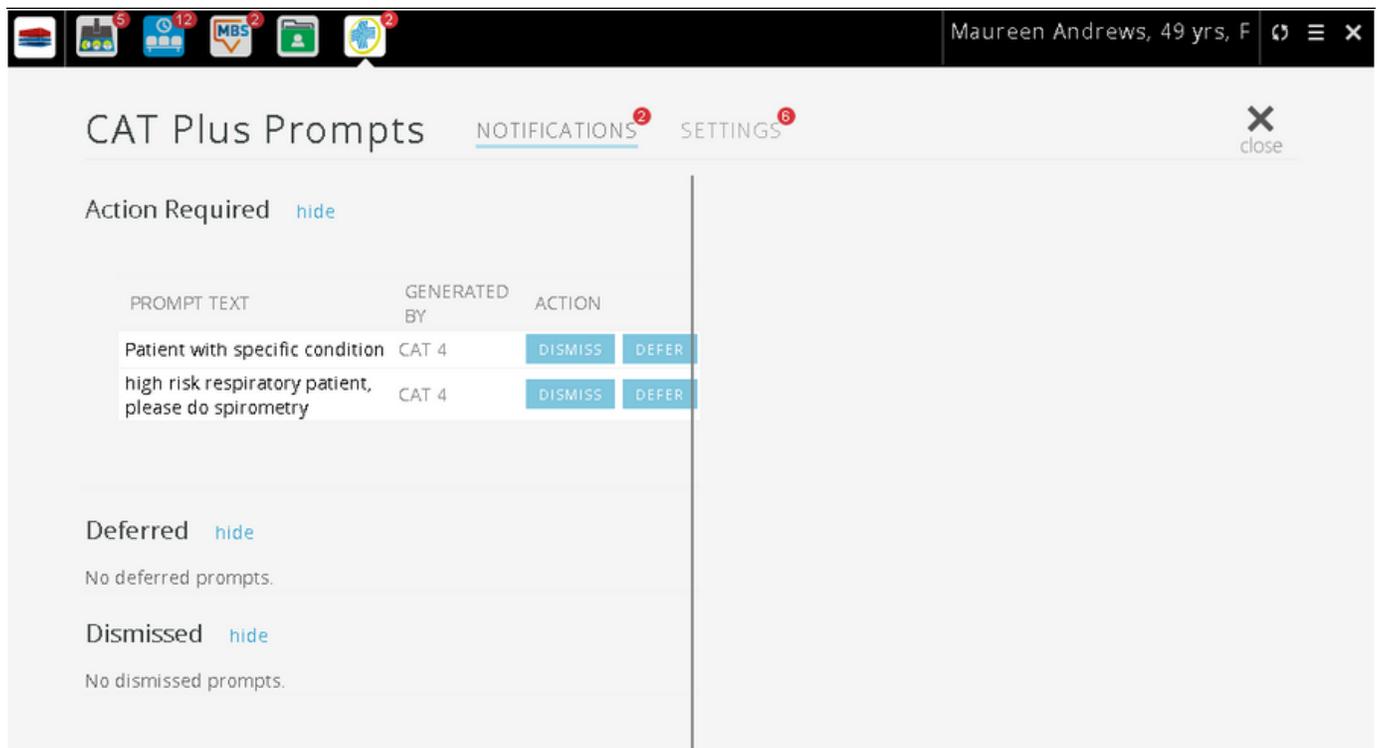

Clicking on the prompt will display further details as shown below. The details show the reason entered when the prompt was created and which reports and/or filters were used. For prompts originating in PAT CAT the organisation who created the prompt is shown along with the date and time. Prompts created in CAT4 will show the CAT4 user who created the prompt. In the 'Reports Used' section only those report sections that match the patient's details are displayed, in the example below the patient is a daily smoker so the 'Daily' is displayed for the 'Smoking' report.

The screenshot shows the 'CAT Plus Prompts' app interface. At the top, there are navigation tabs for 'NOTIFICATIONS' (with a red badge '1') and 'SETTINGS' (with a red badge '6'), and a 'close' button. The main content is divided into two columns. The left column shows a list of prompts under the heading 'Action Required' (with a 'hide' link). The first prompt is 'high risk respiratory patient, please do spirometry', generated by 'CAT 4', with 'DISMISS' and 'DEFER' buttons. Below this, there are sections for 'Deferred' (with a 'hide' link) and 'Dismissed' (with a 'hide' link'), both showing 'No deferred prompts.' and 'No dismissed prompts.' respectively. The right column shows a large blue and green '4' icon. Below it, there are sections for 'Prompt Text' (containing the prompt text), 'Filters Used' (showing 'Has Asthma' with 'Asthma' selected and 'Smoking' with 'Daily Smoker' selected), 'Reports Used' (showing 'Last Spirometry : Nothing Recorded'), and 'Prompt History' (showing 'Creation Date: 5/3/2016 2:26:49 PM' and 'Created By: Matthias').

The prompt(s) will disappear when the patient doesn't meet the criteria specified in the prompt criteria any more. In the example above there are a number of possibilities:


- the patient's smoking status is updated to 'non-smoker'
- the patient's history is updated and Asthma is made 'non-active'
- a spirometry result is entered

Any one of the changes above will make the prompt disappear as the patient now doesn't meet the criteria any more. The same principle works for patients that through a change in their condition or status meet the criteria - these will trigger the prompt to display after their record is updated.


From the cat plus prompts app screen prompts can be dismissed or deferred using the respective button. Dismissing a prompt will open up an optional field for a dismissal reason. Clicking on 'Cancel' will return to the previous screen and not dismiss the prompt.

Dismissed or deferred PAT CAT prompts are recorded and their status will be shown in PAT CAT after the nightly update, for details please refer to the section on [viewing and removing prompts in PAT CAT](#).

Dismissed or deferred CAT4 prompts are recorded and their status will be shown in CAT4 after a new collect, for details please refer to the section on [viewing and editing prompts in CAT4](#).

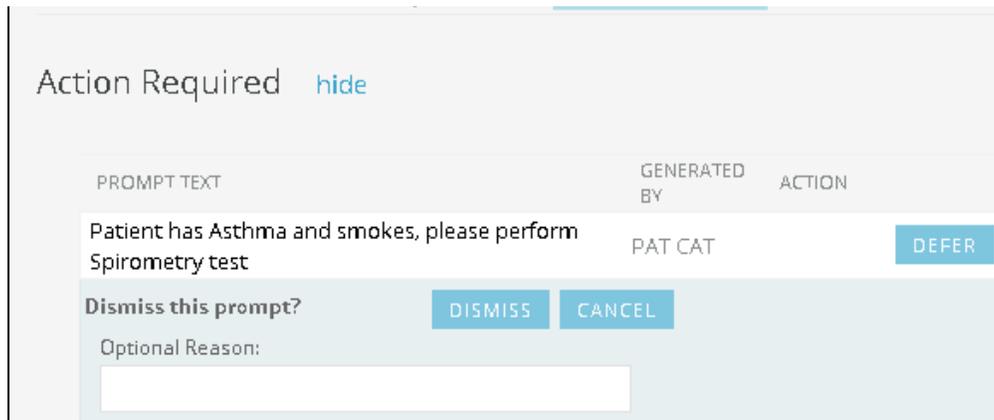

Deferring a prompt will by default defer until three months from today's date, but the date can be changed to a custom date if required. An optional reason for deferring can be entered below.

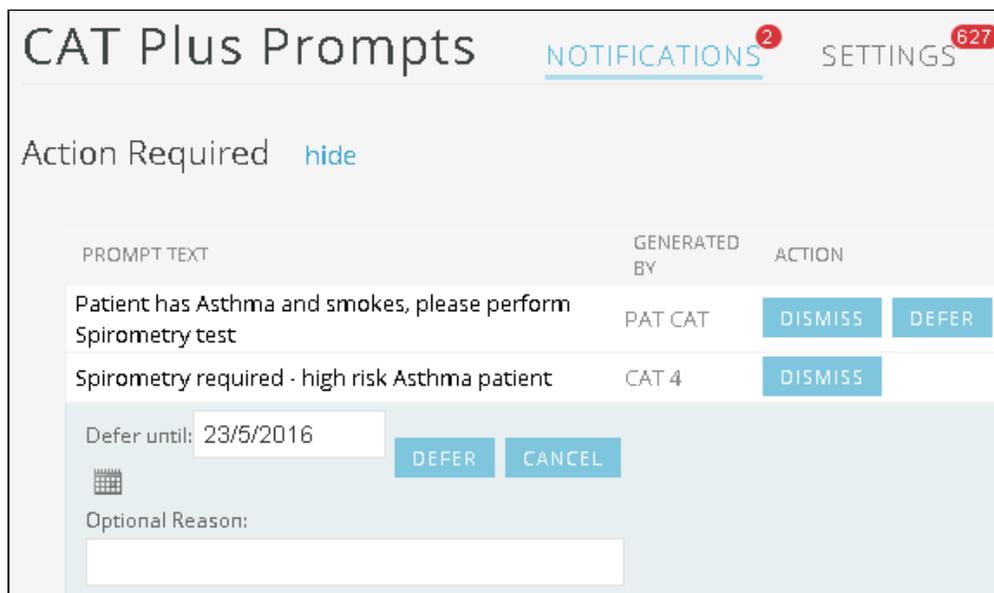

Dismissing or Deferring a prompt will remove this prompt for the currently selected patient only and will affect all users at the practice.

CAT Plus Prompts can be turned on or off depending on the user's preference. This applies to all prompts, regardless if they were created in CAT4 (locally) or PAT CAT prompts created at the PHN or other organisation.

CAT4 prompts that are created at the practice are turned ON for all users by default, where as PAT CAT prompts are turned OFF by default and need to be activated first to be displayed for each user who wants to work with PAT CAT prompts.

To access the settings click on the 'SETTINGS' tab in the cat plus prompts app screen. Existing prompts are displayed sorted by origin, CAT4 prompts first then PAT CAT prompts.

The display status column shows if the prompt is active (green switch) or disabled (grey switch).

